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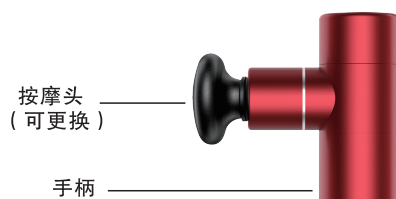
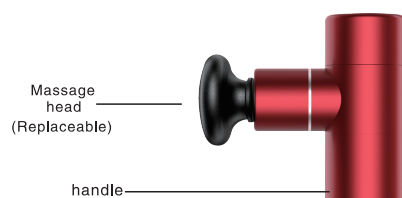
Massage head


Product parameter

-Model: 266
-Rated voltage: 8.4V
-Rated current: 2A
-Input parameters: 5V 2A
-Maximum speed: 3000r/min
-Battery capacity: 1500mAh 2500mHh
-Executive standard: GB 4706.1-2008

Product instructions

On/Off

- Keep pressing the start button for three seconds to turn on or off.
- It will automatically stop after starting 15 minutes, if you want to continue using, please press the button again.

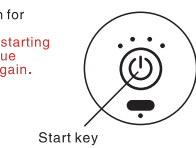
Gear adjustment

- The starting gear is the first gear. The gear is added up one by one when pressing the starting button.
- The highest gear is the fourth one. When in fourth gear, it will come to standby status after you press the button one more time.

NOTE: In order to keep the bone and muscles out of damage, please don't fix point massage for long time. So we set 15 minutes intelligent protection mode.


Replace Massage head

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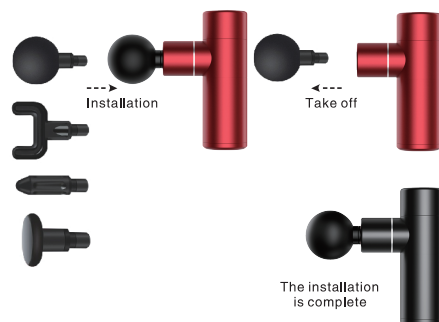
Indicator light, signal light

- The indicator light indicates the size of the fascia gun gear.
- One light for one gear, two lights for two gears, and so on.

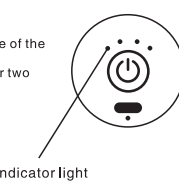
- The signal light flashes once when the machine is turned on.
- The four indicator lights are on once in turn. The red light is always on when charging. The green light is always on after charging.

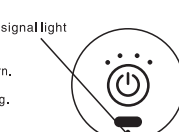- Power on function: the four indicators light up once in turn.
- Power display: below 60%, the signal light flashes twice every 5 seconds. Below 30%, the signal light flashes continuously every 5 seconds. When the voltage is less than 6.5 volts, the warning lamp flashes. If the voltage is lower than 6 volts, it will shut down automatically.

Charging method instructions

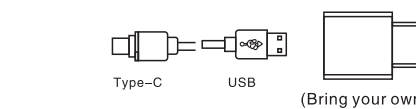- It cannot be turned on during charging. You must unplug the USB charging cable after charging.
- Charge for 2 hours. After being fully charged, it can be used for about 5 hours in the first test. The specific use time will decrease due to the gear and pressure.
- Plug one end of the USB cable into the bottom of the host, and the other end can be connected to the desktop computer. Mobile power, the USB port of the USB power supply is charged. Please note that if you use a USB power supply, do not use a power supply with an output specification exceeding DC 5V. Otherwise, the host will be damaged.

- *To ensure safe use, do not use the fascia gun when charging
- *In order to protect battery life, please charge as soon as possible when the signal light flashes red.

Indication

- Pain and spasm due to muscle contusion, sprain and stain
- To help the swelling flow
- Loose and thickened connective tissue and fascia
- Reduce the accumulation of lactic acid in muscle
- Increase joint mobility
- Eliminate muscle fatigue


Taboo

- Aneurysm, bleeding, blood thinners
- Heart disease, with pacemaker or defibrillator
- Pregnancy, cancer
- The device cannot be used within 3 inches of internal fixation within 90 days after internal fixation
- Sensitive area head, face, cervical vertebra, spine and vertebral body
- The position of tibia, dorsum of foot, dorsum of hand, etc. Near the skeleton with less muscle coverage

There is a lithium battery inside the appliance, so:

- Before discarding the appliance, the battery must be removed from the appliance;
 - When taking out the battery, the appliance must be powered off;
 - The battery should be disposed of safely. Users cannot replace the battery by themselves.
- Take out the battery steps:
-Pry the back cover of the product with a tool and remove it
-Unplug the battery lead and remove the battery.

Checklist

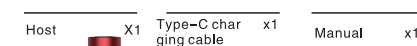Warranty Card Warranty Card

Product Model

Serial Number

Purchase Date

Customer Name

Customer Tel

Customer Email

Warranty service terms

- The product is warranted for one year and the built-in battery is warranted for six months
- The warranty starts from the first purchase date of the product.
- In case of quality problems requiring warranty, please send this warranty card. The purchase documents and products shall be sent back to the dealer or factory for repair.

The following are not covered by the warranty:

- Formal purchase certificate and warranty card with detailed information cannot be provided.
- Damage caused by human or irresistible factors.
- The artificial dismantling of the fascial gun causes irreversible damage.

*The company reserves the right of final interpretation of the above-mentioned after-sale terms and restrictions